
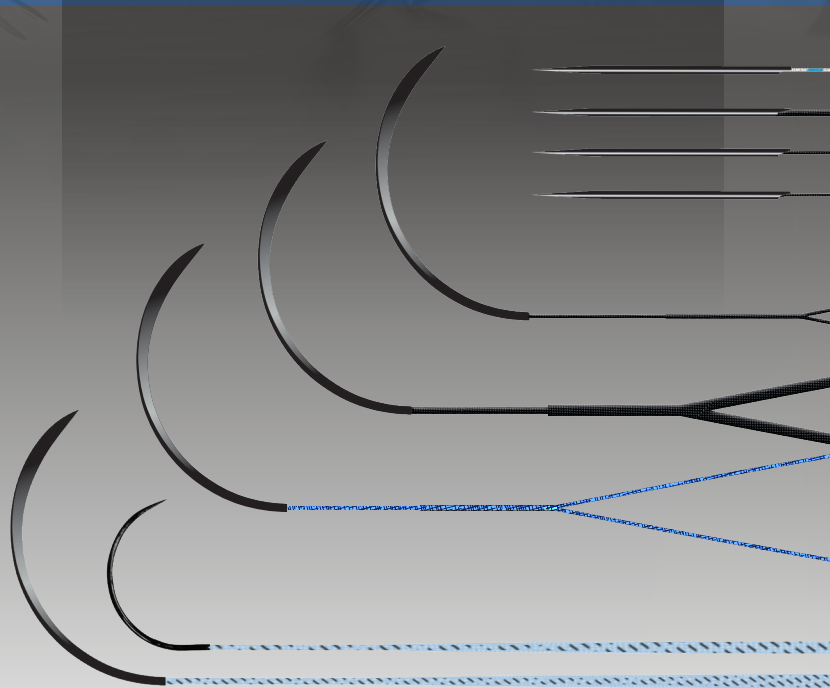


FH LOOP FH LOOP TAPE

FH LINK FH TAPE

NON-ABSORBABLE SUTURES

FH LOOP/FH LOOP TAPE and FH LINK/FH TAPE are non-absorbable sutures composed of a braided ultra-high molecular weight polyethylene (UHMWPE/ Polyester) coated with uniform polybutylene adipate.

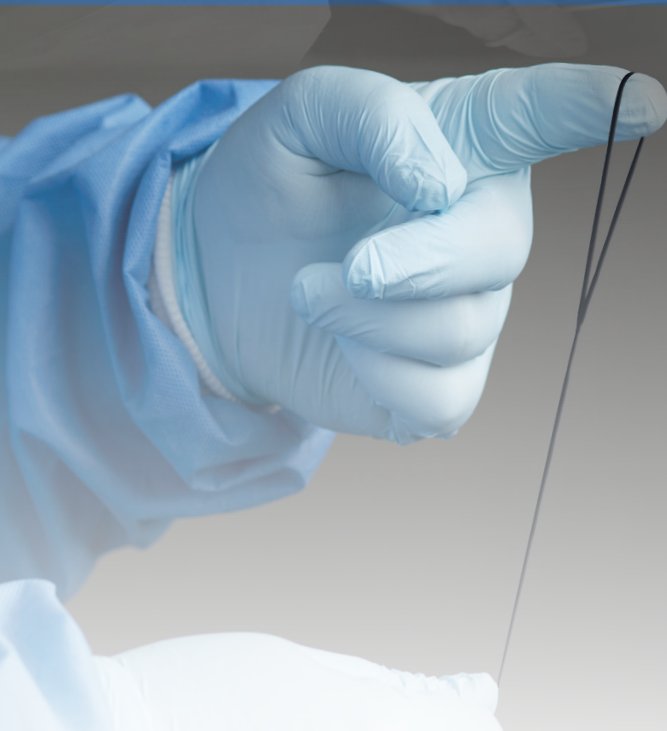

Product characteristics

Needles are composed of 300-Series stainless steel to enhance their penetration force, strength rate and provide good ductility avoiding risks of breakage.

FH LOOP/FH LOOP TAPE and FH LINK/FH TAPE sutures integrate an inter-locking core technology which provides all fiber configurations with a sturdy core within center of the structure to act as backbone for better knot tying and handling characteristics.

Our sutures elicit a minimal inflammatory reaction in tissues, followed by gradual encapsulation of the suture by fibrous connective tissue.

Our sutures are non-absorbable and there are no significant changes in tensile strength retention known to occur in vivo.

Indications


FH LOOP/FH LOOP TAPE and FH LINK/FH TAPE sutures are indicated for use in the fixation of bone and soft tissue in orthopedic procedures requiring ligament or tendon reconstruction.

These sutures are highly appropriate for tendinous graft preparation, ACL repair, rotator cuff repair, Achilles tendon repair, other soft tissue and anchor applications.

Packaging

FH LOOP/FH LOOP TAPE and FH LINK/FH TAPE sutures are supplied sterile, sterilized by ethylene oxide including a black needle for less glare and higher visibility.

Do not re-sterilize.

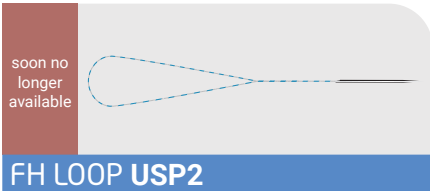


soon no longer available

FH LOOP USP2

- Thread: Ø 0.50 (USP2)
- Length: 60 cm (loop length: 50 cm)
- Colors: blue and white
- 300-series steel (steel and inox alloy) straight, silicone coated and Cobra Black® (black) needle

FH ORTHO reference: 271 197
Supplier reference: HS173

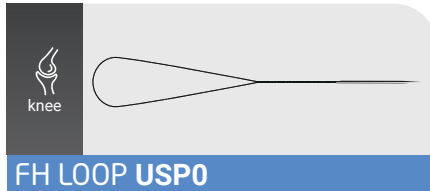


knee

FH LOOP USPO

- Thread: Ø 0.35 (USP0)
- Length: 60 cm (loop length: 50 cm)
- Color: black
- 300-series steel (steel and inox alloy) straight, silicone coated and Cobra Black® (black) needle

FH ORTHO reference: 271 633
Supplier reference: HS345

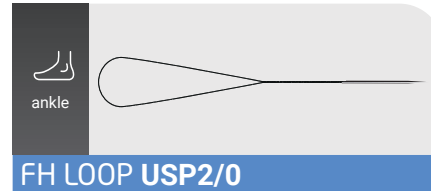


ankle

FH LOOP USP2/0

- Thread: Ø 0.30 (USP2/0)
- Length: 60 cm (loop length: 50 cm)
- Color: black
- 300-series steel (steel and inox alloy) straight, silicone coated, silver needle

FH ORTHO reference: 271 632
Supplier reference: HS343

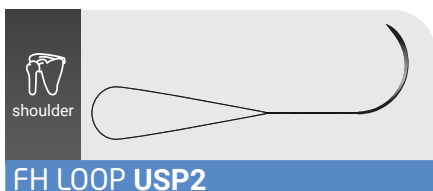


shoulder

FH LOOP USP2

- Thread: Ø 0.50 (USP2)
- Length: 60 cm (loop length: 50 cm)
- Color: black
- 300-series steel (steel and inox alloy), silicone coated, 48 mm black triangular, ½ curved, Cobra Black® (black) needle

FH ORTHO reference: 271 635
Supplier reference: HS179

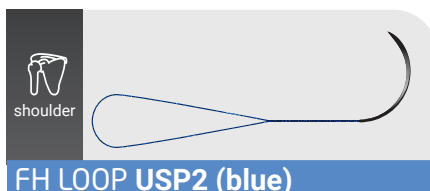


shoulder

FH LOOP USP2 (blue)

- Thread: Ø 0.50 (USP2)
- Length: 60 cm (loop length: 50 cm)
- Color: blue
- 300-series steel (steel and inox alloy), silicone coated, 48 mm black triangular, ½ curved, Cobra Black® (black) needle

FH ORTHO reference: 271 200
Supplier reference: HS266

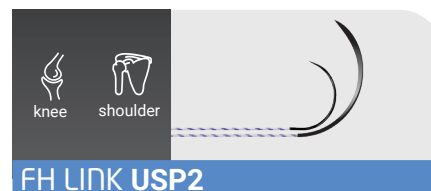

knee shoulder

FH LINK USP2

- Thread: Ø 0.50 (USP2)
- Length: 100 cm
- Colors: blue and white
- 2 Cobra Black® 300-series steel (steel and inox alloy) needles (black), silicone coated:
 - 1 triangular needle, ½ curved, 37 mm length,
 - 1 round needle, ½ curved, 26 mm length.

FH ORTHO reference: 271 203
Supplier reference: HS168




knee ankle

FH LOOP Tape 1.4 mm

- Tape: width 1.4 mm
- Length: 60 cm (loop length: 50 cm)
- Color: black
- 300-series steel (steel and inox alloy) straight, silicone coated and Cobra Black® (black) needle

FH ORTHO reference: 271 631
Supplier reference: HS182

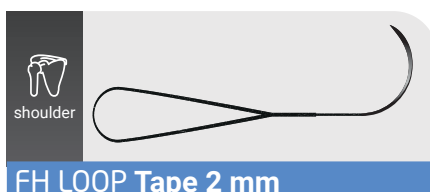

shoulder

FH LOOP Tape 2 mm

- Tape: width 2 mm
- Length: 60 cm (loop length: 50 cm)
- Color: black
- 300-series steel (steel and inox alloy), silicone coated, 48 mm black triangular, ½ curved, Cobra Black® (black) needle

FH ORTHO reference: 271 630
Supplier reference: HS523




shoulder

FH LOOP Tape 2.3 mm

- Tape: width 2.3 mm
- Length: 70 cm (loop length: 60 cm)
- Color: blue
- 300-series steel (steel and inox alloy), silicone coated, 48 mm triangular, ½ curved, Silver

Réf. FH ORTHO : 271 202
Réf. fournisseur : HS519

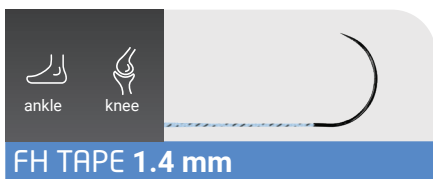


ankle knee

FH TAPE 1.4 mm

- Tape: width 1.4 mm
- Length: 100 cm
- Colors: blue and black
- 300-series steel (steel and inox alloy), silicone coated, round, ½ curved, Cobra Black® (black) needle, 26 mm length

FH ORTHO reference: 271 637
Supplier reference: HS240

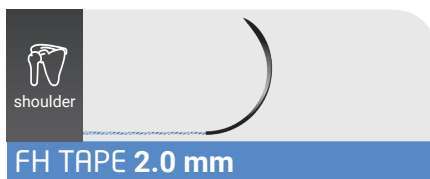

shoulder

FH TAPE 2.0 mm

- Tape: width 2 mm
- Length: 100 cm
- Colors: blue and black
- 300-series steel (steel and inox alloy), silicone coated, 48 mm black triangular, ½ curved, Cobra Black® (black) needle, 37 mm length

FH ORTHO reference: 271 636
Supplier reference: HS524



DISTRIBUTEURS
DISTRIBUTORS

GRUPE
FH ORTHO™

FR, FH ORTHO SAS
3 rue de la Forêt - Zone Industrielle
BP 50009
68990 Heimsbrunn CEDEX - FRANCE
Tél. +33 (0)3 89 81 90 92
Fax : +33 (0)3 89 81 80 11
info@fhortho.com
www.fhortho.com

USA, FH ORTHOPEDICS INC.
OrthoEx
7327 E Tierra Buena Lane
Scottsdale, Arizona 85260 - USA
Phone: +1 (412) 965-0950
customerservice@fhortho-us.com
www.fhortho.com

PL, FH ORTHO POLSKA
Ul. Garbary 95/A6,
61-757 Poznan - POLSKA
Phone: +48 61 863 81 27
Fax: +48 61 863 81 28
biuro@implants24.pl
www.fhortho.com